



# Vitrex<sup>®</sup> Sterilance Press II

## Sterile Safety lancet

**Health workers worldwide face accidental needle-stick injury while caring for patients, putting themselves at risk to blood-borne illnesses including the AIDS virus and Hepatitis.**

**USER SAFETY** needle is fully shielded before and after use.

**USER CONVENIENCE** pre-loaded, pressure activated design.

**PRESSURE ACTIVATION** ensures a good blood sample every time.


**PATIENT COMFORT** high quality, ultra sharp tri-bevel needle or bi-bevel blade, high speed penetration.

**EFFICACY** variety of point configurations and penetration depths to suit most capillary blood requirements.


**APPLICATION** blood glucose and other in vitro diagnostic tests.

**PATIENT SAFETY** sterile, single-use, unable to be re-used.

## Recommended use

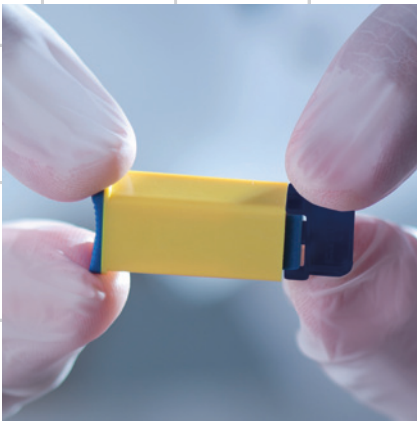

Carefully twist off the protective cap.

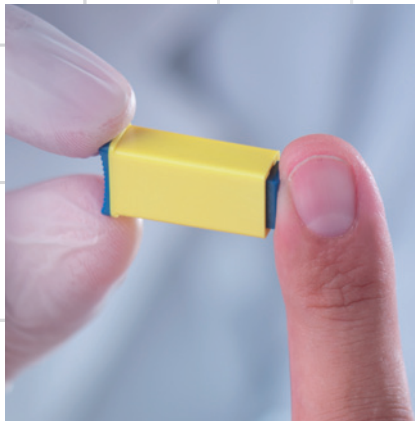

Select site and press the device firmly against the site to release the lancet.



Dispose of safety lancet in a suitable safety container.

**WARNING:** Do not use if protective safety cap has previously been removed.

## Technical Specs.

Colour Code	Configuration	Penetration Depth	Article No.	Application
Green	● 18G blade	1.8 mm	393213	Special baby care and maternity units for heel puncture
Yellow	● 26G needle	1.8 mm	394213	Diabetes clinics/hospital wards for blood glucose monitoring
Orange	● 21G needle	2.2 mm	395213	Haemoglobin (bloodbanks), cholesterol, blood typing, coagulation, blood gas and other diagnostic tests
Pink	● 21G needle	2.8 mm	396213	As Orange but for thicker adult skin
Blue	● 21G needle	1.8 mm	397213	As Yellow but for increased blood flow
Purple	● 28G needle	1.8 mm	398213	New item, less pain

Shelf life 5 years

## Packaging

Primary pack: 100 pcs  
 Transit packaging: 1200 pcs  
 Pallet: 98000 pcs



Transit packaging



Pallet

